

prepFAST

Filtration Station

Automated Sample Filtration System


Compact system that frees laboratory staff from repetitive sample filtration and sample preparation tasks



Features & Benefits

- **Automatic Sample Filtration and Preparation**
 - Eliminates tedious manual sample filtering
 - Functions include filtration, acidification, dilution, internal standard addition and mixing
- **Disposable Sampling Filters**
 - 0.45 micron, 0.22 micron, etc.
 - Automated filter pre-cleaning for low blanks
- **Intelligent Touchscreen Control**
 - User-friendly UI with flexible methods
 - Select your filtration method and click start

Unfiltered Samples Are Made Ready for Analysis

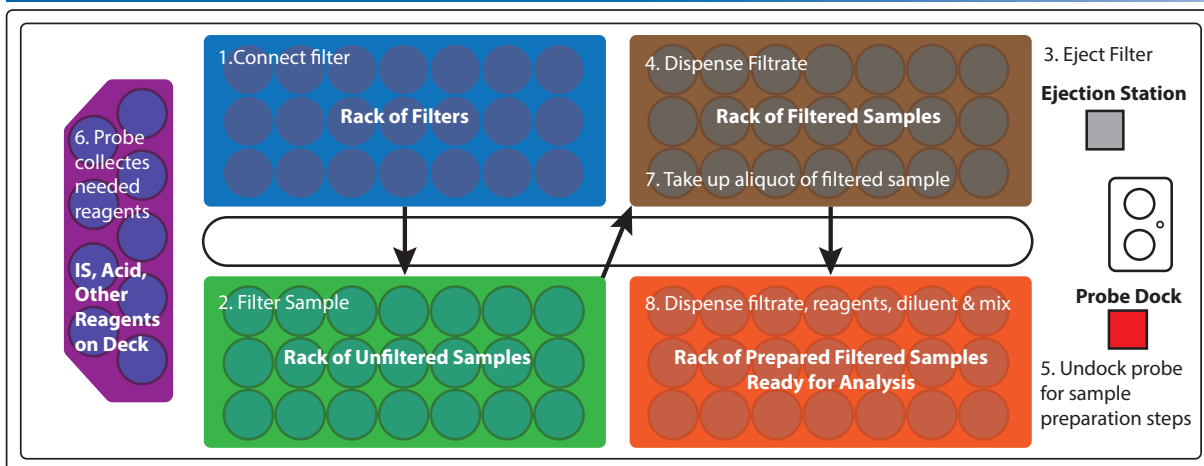

FiltrationStation Features

ULPAClean 10 Filter (Optional)

- Removes 99.9995% of airborne particles with size of 0.1 µm or larger

Clean3V Enclosure (Optional)

- Clean sample environment
- Protects lab from harmful fumes
- Three vertical sliding doors for easy access

AutoAlign Arm

- Self-aligns after encountering capped sample tube

Filter Rack

- 0.45 or 0.22 micron PTFE filters
- Full range of sampling filter options
- 21 or 60 position filter rack options

Reagent Reservoirs

- Accesible 10 L containers (2)

Integrated Touchscreen

- Easy-to-use UI
- Store and access many filtration/prep methods
- Easy customization

Integrated Station Pump

- Rapid and precise dispensing
- Ultra-high purity syringe pump

Dual Rinse Stations

- Ultra-low cross-contamination
- Low rinse consumption

Probe Docking Station



- Demountable sample preparation probe
- Docked during filtration

Filter Ejection Station

- Automatically ejects used filters

Filter Bin

- Collects used filters from the ejection station

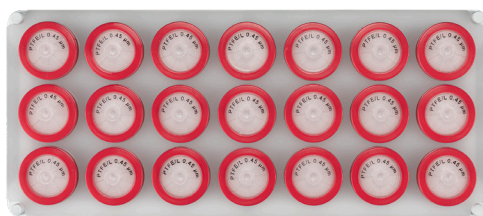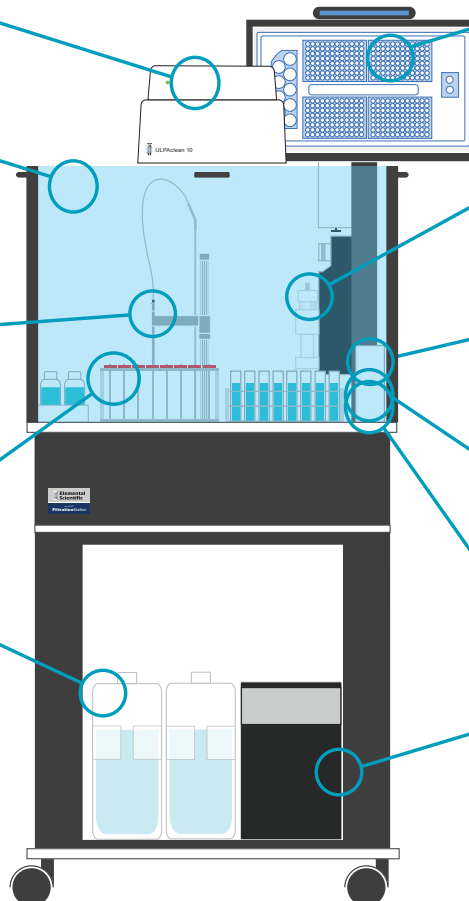

FR21 21 Position Filter Rack



FR60 60 Position Filter Rack

FiltrationStation		
FLTSTN4 4DXCi FiltrationStation		
System Dimensions		
Height	Width	Depth
1753 mm (69")	889 mm (35")	356 mm (14")



© Elemental Scientific | 7277 World Communications Drive | Omaha, NE 68122
 Tel: 1-402-991-7800 | sales@icpms.com | www.icpms.com